



## Coated Allis Bone Clamp – Two Coated Ends

Modification of design by Charles T. Resnick MD

*A traditional Allis Bone Clamp designed with a longer ratchet—for a wider opening to allow a bone and plate to be clamped and locked onto—and coated end(s) to prevent from marring a component surface*


### PRODUCT NO:

1382 [Two Coated Ends]  
Overall Length: 6.125" (15,9 cm)  
Ratcheted Clamp Opens to: 35 mm  
Non-coated-end Width: 4 mm



### PRODUCT NO:

3071  
Overall Length: 8.125" (20,6 cm)  
Jaw Width: .25" (6,4 mm)  
Open Jaw Length: .5" (12,8 mm)

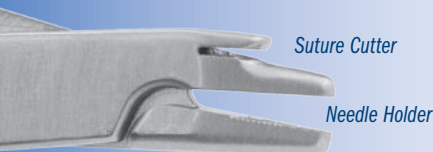

Suture Cutter

Needle Holder

## Bates Needle Holder with Suture Cutter

Designed by James E. Bates, MD



*By trapping the suture and cutting when the forcep is opened, helps to reduce stress on the surgeon's hand*

## Horseshoe Shoulder Frame and Blade Assembly

*Designed to enhance exposure during shoulder arthroplasty procedures*

### PRODUCT NO'S:


2030-00 [Set]  
Includes (1) Frame, (1) Narrow Blade & (1) Wide Blade

### Also available individually:

2030-01 [Horseshoe Frame]  
Overall Dimensions: 7" x 5" (17,8 cm x 12,7 cm)  
Frame Width: .7" (15 mm)

2030-02 [Horseshoe Blade – Wide]  
Blade Width: 22 mm  
Overall Length: 11" (27,9 cm)

2030-03 [Horseshoe Blade – Narrow]  
Blade Width: 14 mm  
Overall Length: 11" (27,9 cm) (Not Shown)

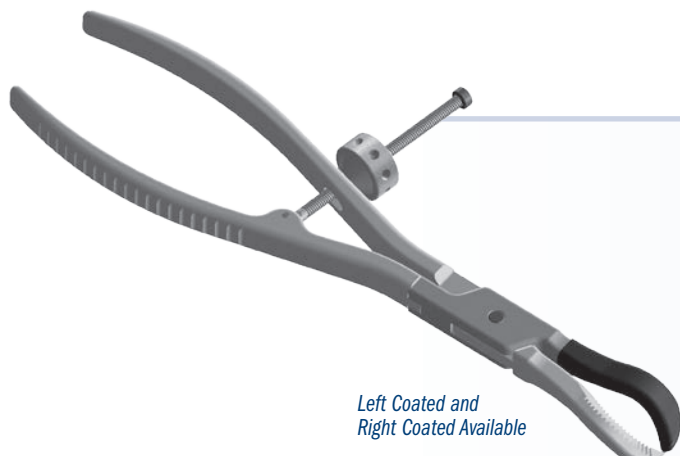

Left Coated and  
Right Coated Available

## Large Bone Clamp with Plate Protection



*Designed to help hold a bone/bone plate in position for reduction—the one-side coated jaw helps to protect from marring the bone plate*

### PRODUCT NO'S:

3659-L [Left]  
Overall Length: 9.125 (23,2 cm)

3659-R [Right]  
Overall Length: 9.125 (23,2 cm)



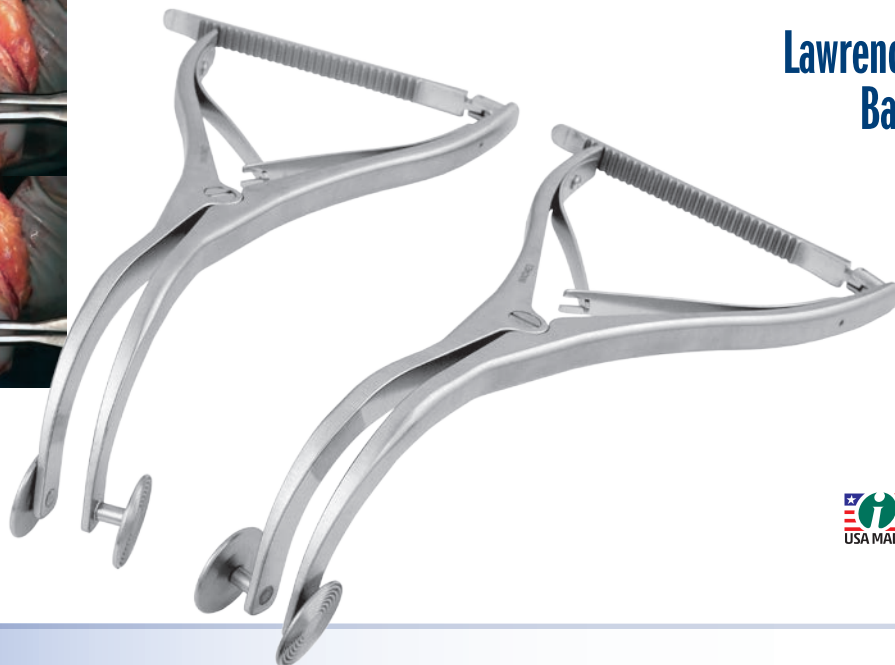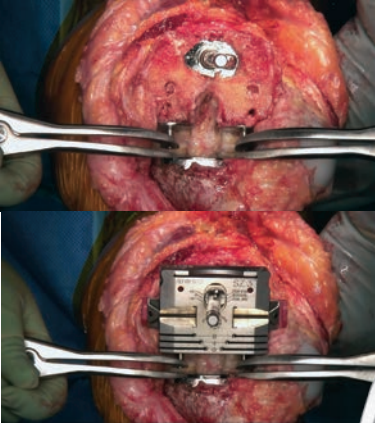

## Lawrence Revision Knee Gap Balancing Tensioner Set

Designed by Jeffrey M. Lawrence, MD

*Designed to help tense the medial and lateral ligaments during total knee surgery*

### PRODUCT NO'S:

1896-01 [Set – Left & Right]

### Also available individually:

1896-01L [Left]

Overall Length: 9.25" (23,5 cm)  
Pad Diameter: 1" (2,5 cm)


1896-01R [Right]

Overall Length: 9.25" (23,5 cm)  
Pad Diameter: 1" (2,5 cm)



## Rudisill Locking Small Bone Reduction Forcep

Designed by Ed Rudisill, MD

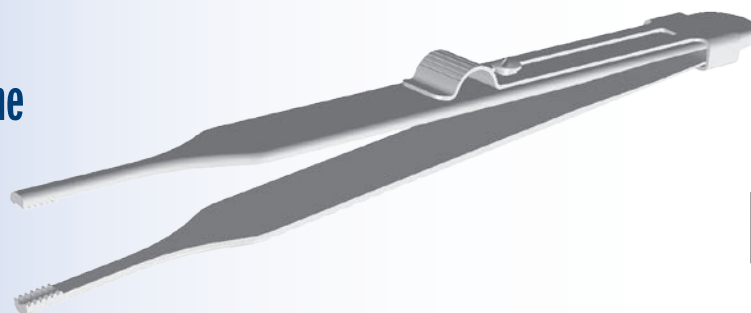

*For reduction of hand phalanx and metacarpal fractures*

### PRODUCT NO:

2017

Overall Length: 6.5" (16,5 cm)



## Stanton Locking Drill Guide Holder for Distal Humeral, Femoral, or Tibial Screws

Designed by John L. Stanton, MD

*Designed to help hold and stabilize a drill guide, allowing the surgeon to obtain 'perfect circles' and drill distal locking screw holes without exposure of the hand to the x-ray beam*

The drill guide is placed over the side of the bone through an incision. The holder is attached to the guide and rested against the skin for stability. With the x-ray on, the guide tube is adjusted by moving the holder until the guide lines up with the hole in the rod. A drill bit is then inserted into the guide.






### PRODUCT NO:

8987

Overall Length: 10.5" (26,7 cm)

Guide Height: 21 mm

Closed Guide Internal Diameter: 5,5 mm

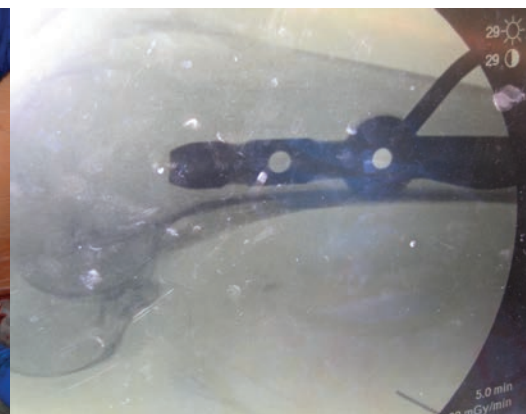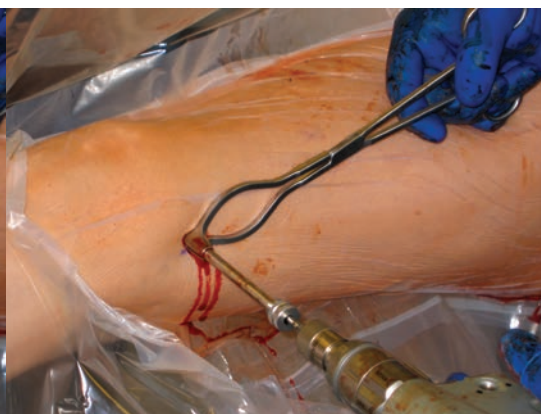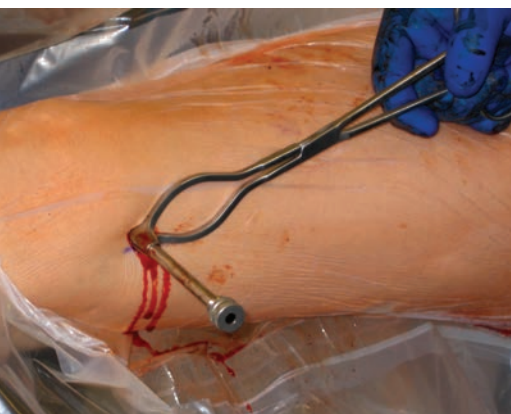

# INNOMED

© 2017 Innomed, Inc.

**Innomed, Inc**  
103 Estus Drive  
Savannah, GA 31404  
USA

Tel 912.236.0000  
Fax 912.236.7766  
[www.innomed.net](http://www.innomed.net)  
info@innomed.net

**Innomed-Europe LLC**  
Alte Steinhauserstr. 19  
CH-6330 Cham  
Switzerland

Tel 0041 (41) 740 67 74  
Fax 0041 (41) 41 740 67 71  
[www.innomed-europe.com](http://www.innomed-europe.com)  
info@innomed-europe.net

**Innomed-Europe GmbH**  
c/o Emons Logistik GmbH  
In Rammelswiesen 9 D-78056  
Villingen-Schwenningen  
Deutschland

Tel 0049 (0) 7720 46110 60  
Fax 0049 (0) 7720 46110 61  
[www.innomed-europe.com](http://www.innomed-europe.com)  
info@innomed-europe.net